


Deep Dive Into My Library

Jamie Clark and Keith Garvey

DIUC22

Making the Deep Dive into My Library

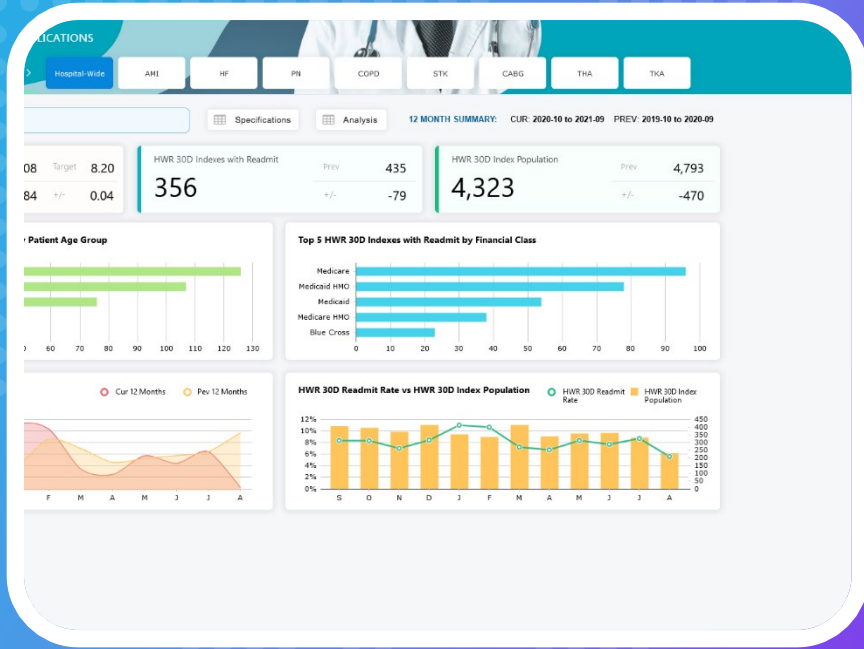


In this session, you will explore My Library and learn how to get the most out of this new tool.

What we will cover today

- ▶ Recap My Library
- ▶ Explore the My Library interface
- ▶ Shortcuts
- ▶ Self Service
- ▶ Visualization Pages (Menu Portlet)
- ▶ Button-style QuickViews

My Library: Interface


Library Bar

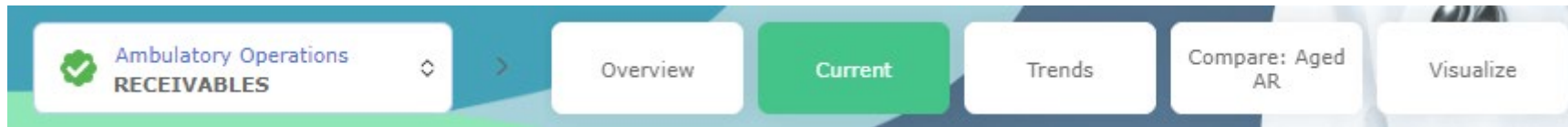

- Show current page and “nearby” displays
- Quick access to all displays through the library button


My Library Popup

- Your whole portal, at a glance
- Sections, shelves, buttons
- Icons for special features
- Easy access to shortcuts and shared folders

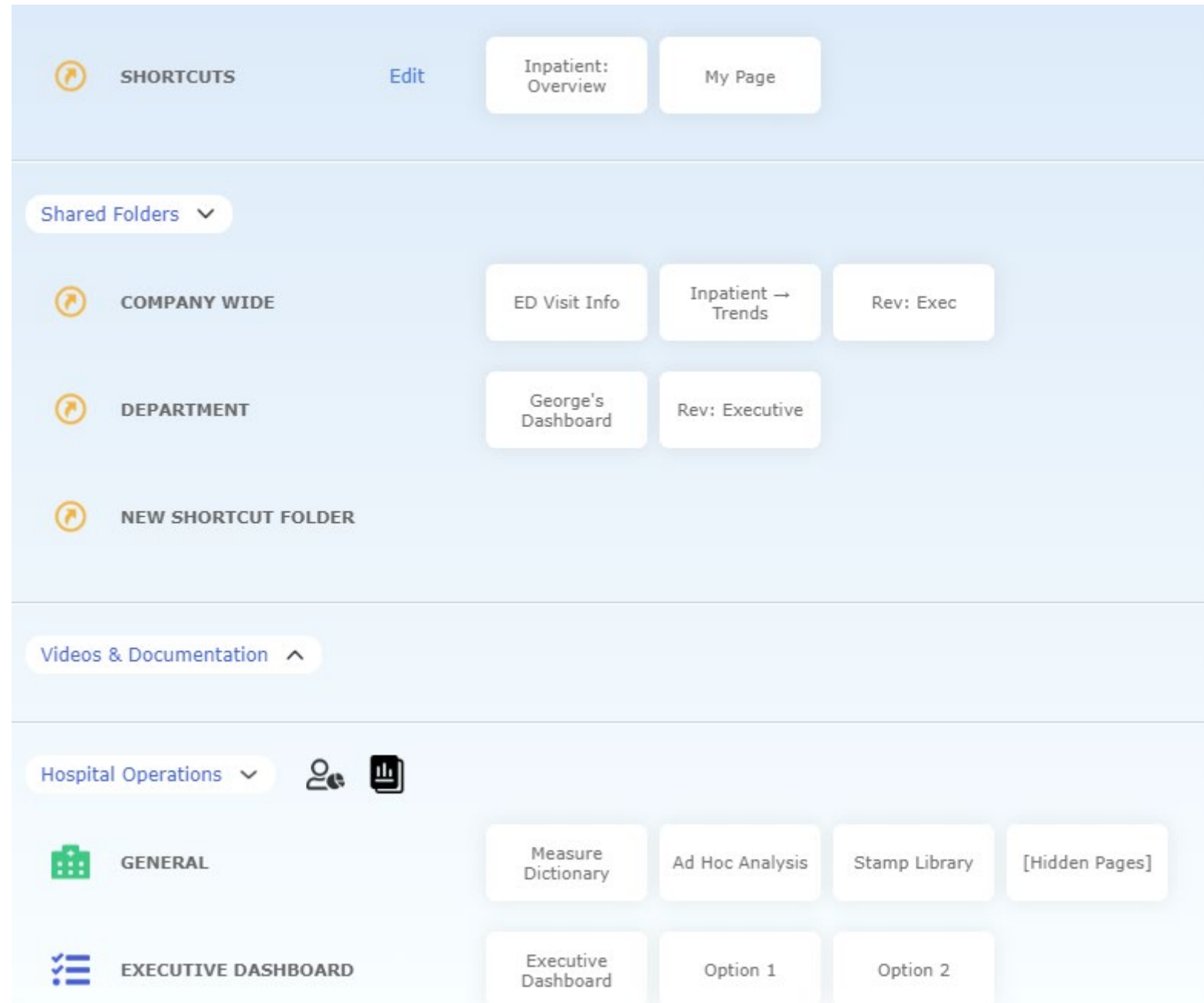

My Library: Menu Portlet "Gallery" mode

- New mode for the Menu Portlet
- Bigger images, better for screenshots
- Special video link option

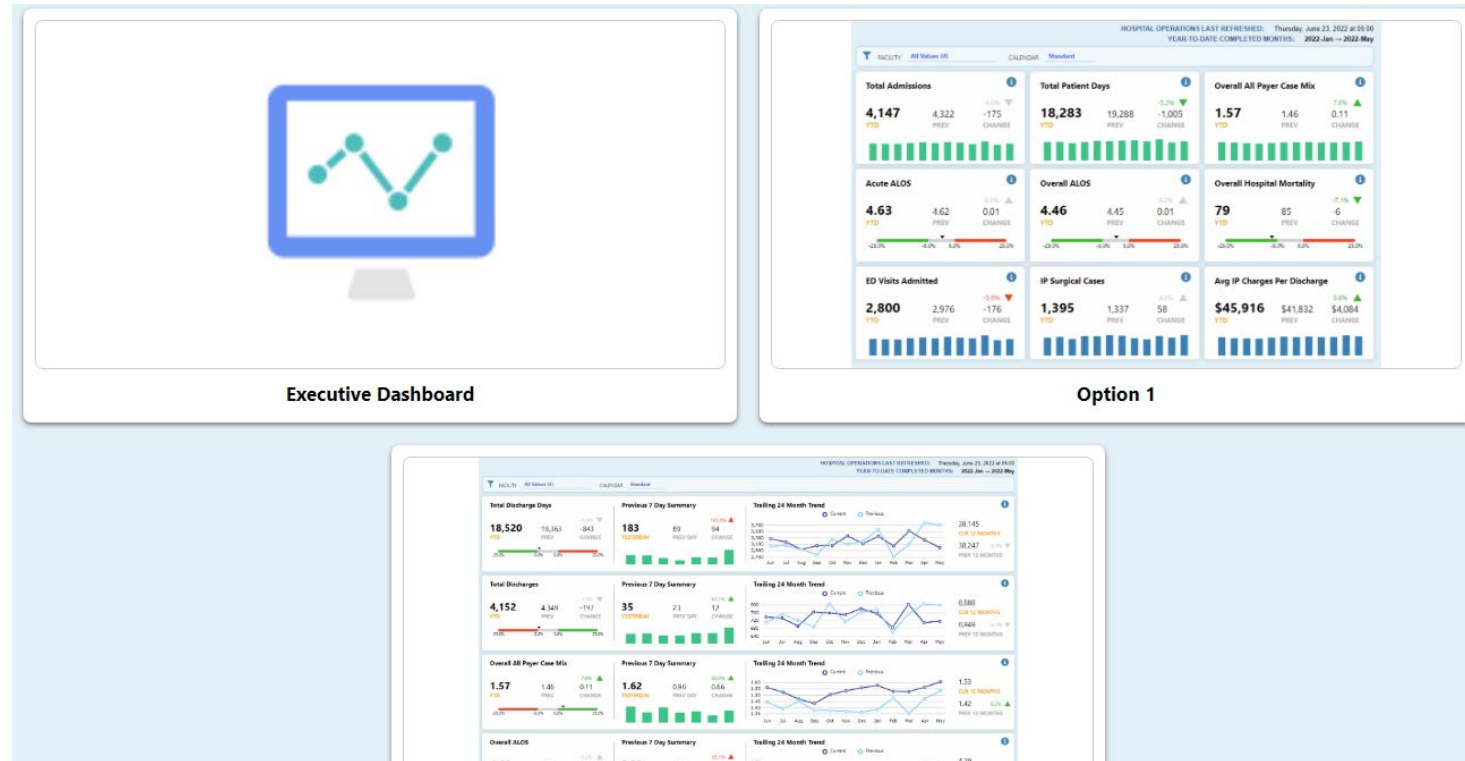


My Library: Shortcuts

- New “Shortcut” system, for use in place of DivePort Bookmarks
- Shortcuts are stored on DiveLine
- Can target a DivePort page or a URL (more to come)
- Easier to reach and manage
- Your first shortcut is your landing page



My Library: Shared Folders

- Create “Shared Folders” on DiveLine
- Groups of users can share shortcuts
- Uses standard DiveLine file access controls

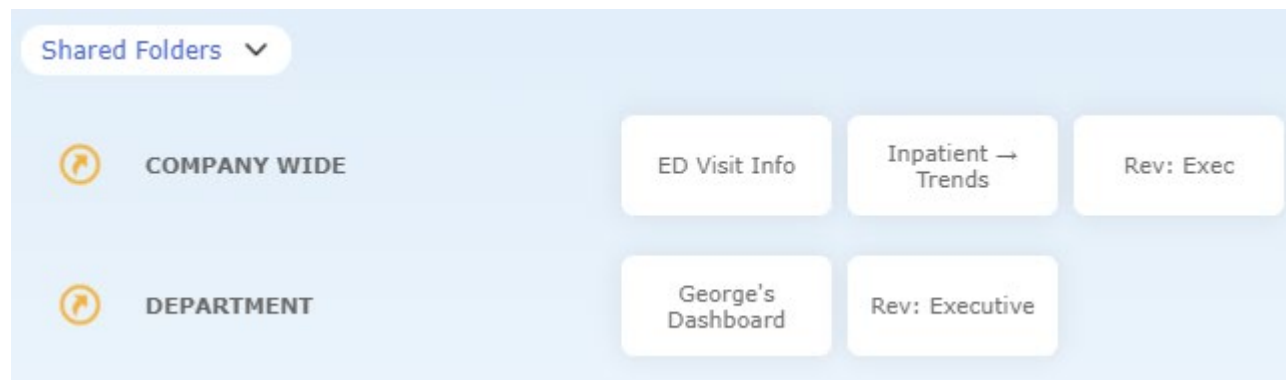

My Library: Simpler Self Service

- New measures dialog in Self Service
- Add multiple measures at once
- New Measure Factory config option: default stamp type
- If set, users can add stamps even faster

Select measures

search

Collapse all

ED Throughput (0/2)

- ☐ ED Admits LOS ED1 Median
- ☐ ED Departs LOS OP18 Median

Hospital Acuity (0/10)

- ☐ Acute All Payer Case Mix
- ☐ Average LACE Acuity Score
- ☐ Average LACE Comorbidity Score
- ☐ Average LACE ED Score
- ☐ Average LACE LOS Score
- ☐ Average LACE Score

ED Volumes (0/6)


- ☐ ED AMA/ Incomplete Treatment
- ☐ ED Visits
- ☐ ED Visits Admitted
- ☐ ED Visits Admitted %
- ☐ ED Visits DC Obs
- ☐ ED Visits Departed

Hospital Population (0/5)

- ☐ AMI Population
- ☐ COPD Population

DivePort: Button-style QuickViews

- Quickly shows all options to the user
- Customizable via skin properties
- Can allow the user to modify the display with the click of a button

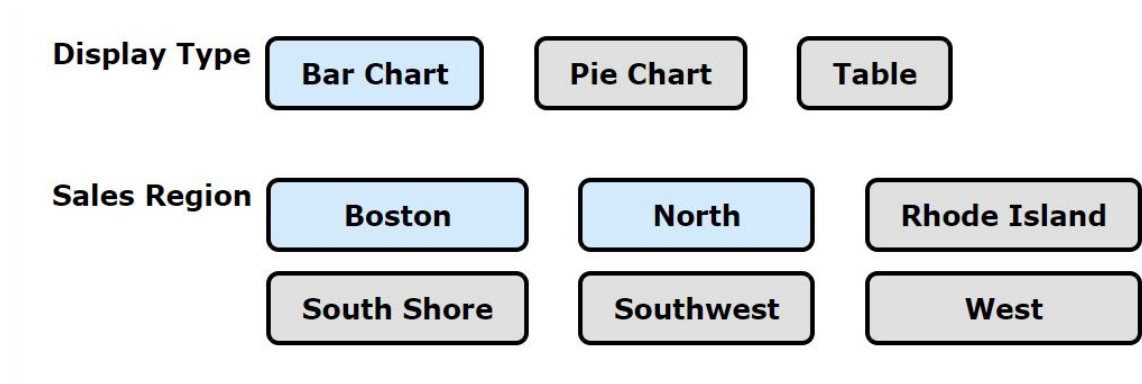



My Library Details

Turning it on

- Able to upgrade to 7.2 without My Library
- Enable in Portal Settings as the default interface
- Or set on a per-environment basis for testing
 - Avoid mixing interfaces in a live site





My Library Details

Sections

- Sections in the popup are Environments
- Start collapsed or expanded
 - Customize in environment settings
- Your current session's collapse state is remembered
- Resets when you log off



My Library Details

Inside Sections: Shelves, Icons, Buttons

- Each "shelf" is a top-level page
- Sections can have icons leading to special pages
- Each shelf has an icon
 - Customize in page settings
- Each shelf button is a second-level page





My Library Details

Administration

- Use the sidebar to create and move pages
- Must reorganize to apply My Library to an existing site



Shortcut Details

Bookmarks vs Shortcuts

- Bookmarks are stored in DivePort's portal database
- Shortcuts are stored in files in DiveLine's Home Project
 - Will later be used in products like DiveTab and ProDiver
- Bookmarks only point to DivePort pages
- Shortcuts point to Markers, DiveTab pages, external URLs, etc.





Shortcut Details

Migration

- Enabled in portal settings
- DivePort automatically migrates your bookmarks to shortcuts

Requirements:

- DiveLine Home Project
- Home Directories for each user



Shortcut Details

Creating Shared Folders

- Make a folder in a special place in the Home Project
- Set access controls using DiveLine access rules:
 - If a user can read the folder, they will see the shortcuts
 - If a user can write to the folder, they can create and edit shortcuts there





Visualization Page Details

New Menu Portlet feature:

- “Gallery” mode
- More space for each button → better for screenshots



Visualization Page Details

What is a “Visualization Page”?

- Contains a Menu Portlet in gallery mode
- Contains screenshots of highlighted portal pages
 - Originally meant to be “live” but security is a priority
 - Requires generating an image for every page for every user
- Another way to help users find useful content
- Optionally, add helpful video tutorials





Button-Style QuickView Details

New DivePort QuickView display type

- Useful for QVs with a fixed or small number of values
- Fewer clicks to change values
- Single and multiple-selection modes





Button-Style QuickView Details

- Works for Dimension QuickViews
 - Be aware of dimensions with many values
 - More values might affect page layout
- Even more effective for List QuickViews